

ZipChip

Identification of adeno-associated virus capsid proteins using ZipChip CE/MS


Yun Zhang, Yan Wang, Zoran Susic, Li Zang, Svetlana Bergelson, Wei Zhang

Anal. Chem., Just Accepted Manuscript

DOI: doi.org/10.1016/j.ab.2018.06.006 Publication Date (Web): June 8, 2018

Copyright © 2018 American Chemical Society

ABSTRACT: A simple and rapid identity test of adeno-associated virus (AAV) serotypes is important for supporting the AAV gene therapy development, as it relates to its efficacy and safety. The current mass spectrometry-based identity tests require extensive sample preparation steps, relatively large sample quantities and long analysis time. Herein, we describe a simple and novel microfluidic ZipChip CE/MS method used to characterize AAV capsid proteins. The three capsid proteins of AAV2 were separated and identified within 4 min using 5 nL of sample directly from a polysorbate-containing formulation buffer. This rapid method can be suitable to confirm AAV serotype identity.

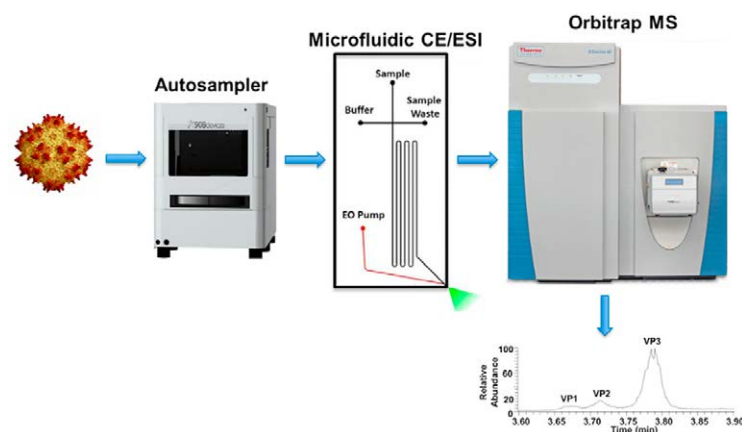

[Click here](#) to view the entire article.



 **908devices**

www.908devices.com

ZIPCHIP IS FOR RESEARCH USE ONLY

Zip Chip is subject to export controls including those of the Export Administration Regulations of the U.S. Department of Commerce, which may restrict or require licenses for the export of product from the United States and their re-export to and from other countries. Patented technology www.908devices.com/patents © 2019 908 Devices